

impress yourself.



Skinfresh 6



GB Operating instruction



IP20



A946

Druck Nr. 29494014en / 10.12

English

Index

Interesting information on collagen	3
Important notes	3
Information on UV exposure	4
Initialization	5
After exposure	5
Memo disc	5
Period of use of the lamps	5
Lamp change	6
Repair	6
Technical data	6
Service address	7
Guarantee	7
Handling procedures for return shipments (RMA) - Details for all returns !	8

Dear customer,

With the purchase of a DR. KERN collagen unit, you have selected high quality technology and a high-performance unit. During production of this unit, we have taken all measures possible to guarantee fault-free, safe operations.

Furthermore, you can make a significant contribution to ensuring that you are satisfied with your unit for a long time. Please read the operating instructions carefully and assemble the unit following the enclosed assembly instructions. The full function of the unit can only be guaranteed if it is assembled correctly.

Interesting information on collagen

The collagen exposure uses visible red light with a frequency of 633 nm in order to stimulate the body's own process of "Photo Bio Stimulation". This process is mainly responsible for the development of collagen and elastin and provides a more taut, youthful and fresh complexion. The penetrability of the skin for cosmetic products is also increased so that the effective nourishment can penetrate deeper into the skin. **Skinfresh** can be used on the whole body and is particularly suitable for the: face, eye contours, hands, feet, thighs and legs.

- Regeneration and rejuvenation of the skin
- Visible anti-aging effect: combats wrinkles and fine lines.
- Improves the skin's appearance for uneven skin
- Clinically-proven stimulation of collagen and elastin
- Preparation of the skin for cosmetic applications
- Reduces the effects of damaging environmental influences

Important notes

 **DANGER**



Read the instructions carefully.

When using a collagen unit, it is necessary that possible overreactions of the skin are carefully observed and the following advice is precisely followed.

Please consult a physician before use if the below points apply to you.

- Persons who suffer from sunburn;
- Persons with more than 20 moles on their bodies;
- Persons who have unusual moles - consult a physician in case of any doubt;
- Persons with unusual discolored areas of skin;
- Persons who suffer, or have suffered from skin cancer or persons who are predisposed to this;
- Persons whose close relatives have been diagnosed with skin cancer;

DANGER

- This unit is not intended for use by persons (including children) with restricted physical, sensory or mental capacity or by persons without the necessary experience and/or knowledge unless they are supervised by someone responsible for their safety or receive instructions from them on how to use the unit.
- Children do not recognize the dangers that can occur when using electrical equipment. Therefore, never leave children to work with electrical equipment when unsupervised.
- Check whether the voltage information on the unit label matches the mains voltage in your building before connecting the unit to the mains.
- Do not use the unit in bathrooms, humid areas or near to swimming pools.
- Do not switch the unit on in areas at risk of explosion.

The collagen unit should be connected to the electricity mains in accordance with the technical data .

CAUTION



Air intake slits or hot air ventilation openings may not be covered, otherwise the unit will overheat.

In new units, the lamps may flicker or display shadow effects. This does not have an effect on the function.

Slight blackening may appear on the ends of the lamp but this does not have any influence on the function of the unit or the life span of the lamps.

CAUTION

Collagen lamps may only be replaced with types marked on the unit.

Repairs on the collagen unit may only be carried out by the factory customer services.

CAUTION



The collagen unit should be disconnected from the mains during any work (lamp replacement or cleaning work) on the combination unit (lamps, starter replacement or cleaning work).

Information on UV exposure

- Please observe the recommendations with respect to exposure periods, exposure intervals and distance from the lamps.

This is why the following cautionary measures must be taken:

CAUTION

- Remove cosmetics before exposure.
- Seek medical advice if stubborn swellings or wounds form on the skin or if pigmented moles change their appearance;
- Protect sensitive areas of skin such as scars, tattoos and sexual organs before exposure;
- Please consult a physician before applying further exposure if side effects, such as itching, occur within 48 hours after the first use.

Initialization

Place the unit on a stable surface (table) and connect to 230 V alternating current.

Set the exposure time with the rotating button on the time switch; the unit will be switched on.

After setting the time, make yourself comfortable keeping your face at least 30 cm away from the unit.

Once you have set the exposure time, the exposure begins.


After expiry of the set exposure time, the unit automatically switches off.

After use, pull the plug from the socket.

Before putting the unit away (after use) allow it cool for 10 minutes.

CAUTION

Maintain a minimum distance of 30 cm from the face.

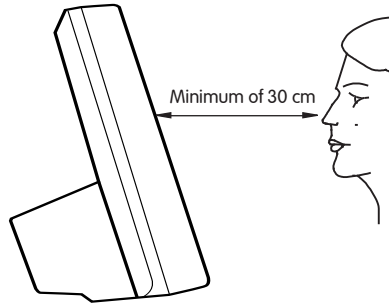


The first exposure should be in accordance with your own level of comfort. The skin should not feel tight or react allergic or over-sensitively afterwards.

Recommended use: 2-3 time a week. Each approximately 20 Minutes.

After exposure

Use a moisturizing lotion after the session to care for your skin. This will keep the appearance of your skin fresh.

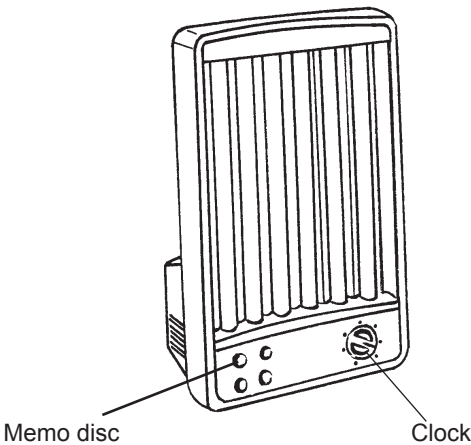

Observe the information and follow the exposure instructions.

Memo disc

So that you know how frequently you have used the unit, move the memo disc along by one number after each use. The other three memos are intended for other people.

Period of use of the lamps

The collagen fluorescent lamps have a maximum period of use of 500 hours. After this time, the lamps still illuminate but the collagen exposure is reduced so that the lamps should be replaced. Blackening of the lamp ends does not affect the function and the service life of the lamps.

Lamp change

Please only use original collagen lamps.

CAUTION



It is compulsory to disconnect the mains plug from the socket.

- Loosen the four fixing screws on the screen on the back of the unit and take the screen off.
- Rotate each lamp so far that lamp holder slits can be taken out from the front.
- Insert new lamps by sliding the contact pins up and simultaneously sliding down into the lamp holder slits, then turn the lamps a quarter turn on their own axis.

Repair

CAUTION



It is compulsory to disconnect the mains plug from the socket.

- Repair work on and in the unit should only be carried out by an expert.

Our service address can be found on the same page as the guarantee terms.

Technical data

Mains voltage: 230 V
Lamp power: 25 watt
Lamp length: 44 cm
Number of lamps: 6
Power consumption: 160 watt
Weight approx.: 5 kg

Technical changes are reserved

Guarantee

The guarantee is taken over according to the legal regulations at present.

Manufacturer's warranty

- The period of warranty starts from the date of purchase and is in principle valid for 12 months.
- Always include the completed warranty certificate when returning equipment.
- The warranty expires for appliances which have been modified without manufacturer's explicit agreement.
- Damages caused by incorrect operation or handling through non-authorized persons are not covered under the terms of warranty.
- In the event of a claim, please indicate the serial number as well as the article code number and type name with expressive description of the fault.
- This warranty covers damaged parts but no defects due to wear and tear.

In case of complaint please return the equipment in its original packaging or other suitable packaging (caution: danger of transport damage) to our service department.

Always include the completed warranty certificate when returning equipment.

Possible shipping costs arising from the transport to and from point of repair cannot be borne by us.

Outside of Germany please contact your specialist dealer in case of warranty claims. Direct warranty processing with our service department is in this case not possible.

Equipment start-up date:

Stamp and signature of the authorized electrician:

Service address

K. Oppermann GmbH
Wilhelmshäuser Straße 81
D-34346 Hann. Münden
Germany

Tel.: +49 (0)5541 95563-20
Fax: +49 (0)5541 95563-22

E-Mail: info@oppermann-online.de
www.oppermann-online.de

Please keep this address in a safe place together with the installation guide.

To help us answer your questions quickly and competently, please provide data printed on the ID plate, to include system type, article no. and serial no., in all inquiries.

Handling procedures for return shipments (RMA) - Details for all returns !

Dear customer

we hope that you will rejoice in the ordered articles. Just in case that you are not entirely contented as an exception, please follow the procedures specified below. This enabling us to ensure a quick and smooth handling of the return shipment.

Please absolutely respect for all returns!

- Please add the available **RMA-voucher** always **completely filled out** together with an **invoice copy** to the return shipment! Do not stick it on the goods or on the packaging. **We do not accept the return shipment without these papers.**
- Not prepaid parcels **will be refused** and returned to Sender! Please always ask for the **RMA-No.** for the cheapest return.
- **Please pay attention that** the goods have to be sent back **without visible marks of use** in the **original scope of delivery and in original packing.**
- We recommend to use an **additional solid and break-proof covering box** which should be padded out with styrofoam, paper or similar. Transport damages as a result of faulty packing are for the sender's account.

Form of complaint:

1) Transport damage

- Please check the content of your parcel immediately and advise the forwarding company of a claim (parcel service/freight forwarder)
- Do not use damaged goods!
- Ask the forwarder **for a written acknowledgement of the damages.**
- **Report the claim promptly by phone to your dealer.** He will discuss with you how to act in this case.
- If the transport box has been damaged, please use an additional covering box. Do not forget to add the **acknowledgement of the damage of the forwarding company !**

2) Faulty goods

- The implied warranty period is 2 years. Please contact your dealer in case of **faulty or wrong articles or missing accessories.** He will discuss with you the individual case and try for immediate and customer-friendly solution.
- For **economic returns** within Germany you will get an **RMA-number from the manufacturer.**
- All returns have to be in the **original packing of the goods with corresponding accessories.** Please repack the goods to avoid damages. In case of wrong delivery, please do not use this article !

3) Problems of installation and functioning

- Please **read the manual carefully first of all** and pay attention to the indicated assembly or installing instructions.
- **Your dealer should be the first contact person** because he knows his products best and also knows possible problems.
- **In case of function problems with an article,** please check at first whether there is an obvious material defect. The quality system in our factory reduces malfunctions of new appliances to almost zero.